

LEGS BEAUTICIAN

LEGS BEAUTICIAN



CE

OPERATION INSTRUCTIONS



CONTENTS

- 01 Function of Product
- 02 Using Effect of Product
- 03 Part name of Product
- 04 Operation Instruction of Computer Panel
- 05 Using Method
- 06 Safe attention ions
- 07 Product Nursing

First , we appreciate your support .Please read the instructions carefully before operating this machine , and well keep it for later use . Note: we reserve the right to revise or explain the design of this product without further notice.

LEG BEAUTICIAN SPECIFICATIONS

- **MODEL:LF-288 Series**
 - **VOLTAGE:AC220-240V~50Hz 1.6A**
 - **POWER:80W**
 - **RATED WORKING TIME: 15 MINUTES**
- A/C ADAPTOR
- **INPUT:AC220-240V~50Hz 1.6A**

SAFETY

○ Notes about power plug and cord.

Please check whether the voltage of power supply is fit adapter standard, Do not disconnect the AC plug from the electrical outlet roughly. DO not connect or disconnect the AC plug from the electrical outlet with wet hands . It may result in electrical shock or machine trouble . When disconnecting the AC plug ,pull it not by holding the plug and not the cord. Do not place the AC cord under the massager or any heave object. Do not wind the AC cord around the massager as it may damage the cord and result in fire or elecnal shock . Do not operate the massager with a damaged cord or tamper with the cord, contact the service centre for inspection. Please use the adapter and attachment we recommend . Do not use the massager if the electrical outlet is excessively loose.

○ Product care and storage.

Be sure to remove the AC plug from the electrical outlet before cleaning the massager. Do not use benzene ,thinner or any other abrasive solvent on your massager. Use a dry cloth when cleaning the massager.Never use a damp cloth to clean the massager or the power switch. When necessary ,detach the fabric covers for cleaning to maintain good hygiene. When not in use for a long period of time, it is advisable that the AC power cord of the massager be coiled and the massager stored in a dust and moisture free environment. Do not store the massager near intense heat or an open flame .Do not leave the massager under direct sunlight for long period of time.

SAFETY

○ Usage Environment

Do not use the massager in an excessively humid environment such as the bathroom.
Do not use the massager on a wet body, this may result in a short circuit or an electric shock.

○ Suitability of User

The following users must consult their doctor before use
Those suffering from cancer.
Those suffering from osteoporosis.
Those suffering from malignant tumors
Those who have a high fever
Women who are pregnant or menstruating

○ Safety Precautions during Use

Do not use the massager for any purpose, other than massage, this may result in accident or injury
Do not use the massager when the fabric cover is torn, or parts of it are missing.
Do not let children, or persons who are not capable of handling equipment, operate the massager. please keep the massager out of reach of children.
Do not stand or sit on the massager.
Do not insert metallic objects, or flammable materials into the main body. this may result in electrocution, short circuit and possible fire.
Do not use the massager within an hour of mealtime.
Do not cover the massager during operation. This may result in a rise in motor temperature and subsequent breakdown or fire.
Do not unplug switch off the massager during operation. This may result in injury.
Do not sleep while using the massager.

○ Instances when the massager should not be used

In the following events, please switch off the power immediately and disconnect the AC power cord from the electrical outlet

If water is accidentally spilled on the massager, this may result in electrical shock.
If severe or abnormal pain is felt during massage. stop operation immediately and consult your doctor.
If you detect any malfunction or any other abnormal condition during operation.
If there is a power failure. injury may occur when power is restored unexpectedly.
If there is lightning

○ Assembly and repair the massager.

The massager must only be serviced or repaired by the service centre. users must not disassemble or repair the massager themselves.


FEATURES

1. An entirely new massager specifically designed to massage your feet, ankles and calves all at the same time.
2. 4 parts of kneading discs in one appliance to massage your feet, ankles and calves.
3. Flexible kneading discs create a professional massage experience.
4. Strong kneading action penetrates deeply into the muscles of your feet, ankles and calves.
5. Pressure nodes located under your feet simulate the reflex points on the soles of your feet.
6. Personal-preference of 3 kneading massage programs and 3 vibration modes.
7. The synergistic combination of kneading and vibratory action gives you a most relaxing and revitalizing massage experience, while achieving a beneficial reflexology and beautifying effect.
8. Strategically located buttons on the switch panel allow you to control the massage function easily.
9. Auto shut-off feature activates after the massager has been in use for 15 minutes.
10. A pair of detachable fabric cover (with zip) to allow for cleaning to maintain good hygiene.
11. **DO NOT accelerate circulation and USE THE HEATED FUNCTION ACCORDING YOUR NEED.**

EFFECTIVENESS


1. The kneading and vibratory actions go deep into the different layers of muscles of your feet, ankles, and calves, leaving your rejuvenation and relaxation.
2. Kneading and vibration massage target vital reflex points on the feet and soles, achieving positive reflexology effects.
3. Stimulates blood circulation, thereby preparing and warming the muscles.
4. Relieves tension and reduces stress.
5. Relieves sore and tired feet effectively.
6. Helps to tone, shape and beautify your legs.

NAMES OF PARTS AND FUNCTIONS

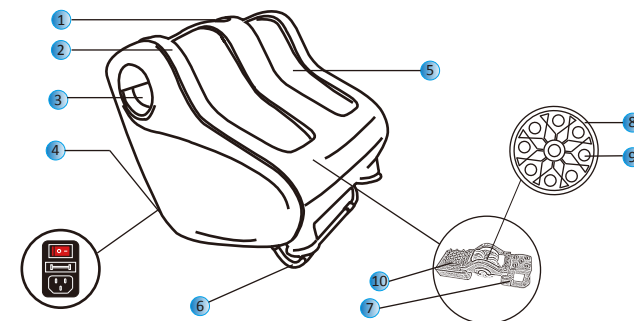



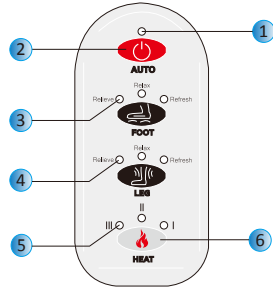

1. Computer panel
2. Shell of mainframe
3. Handle
4. Power switch (Back)
5. Dot cover (with zipper)
6. Bracket
7. Roll subassembly
8. Roll ball
9. Magnet
10. Heat up part

Operation Instruction of Computer Panel




1. indication light
2. power button
3. indication light of foot massager speed
4. indication light of leg massager speed
5. indication light of Heat
6. Heat button

To activate the kneading discs only


1. Place the lower section of your legs into the fabric-covered slots. See picture .
2. Press the  button and the default mode will be "Refresh", press it again or twice, and mode will switch to "Relax" or "Relieve". Different kneading modes for each level.
3. Press the  button and the default mode will be "Refresh", press it again or twice, and mode will switch to "Relax" or "Relieve". Different kneading modes for each level,
4. Press the  button to start function of heating
5. Press the  button to stop the kneading massage.

Cue


Period of any function is 15 minute,time again when change function.

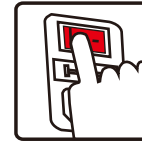


Press the power switch in the back of massager, the indicating light of stand-by will turn on.

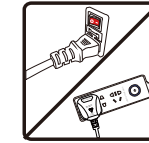


insert the adapter to the massagers back socket by correct direction, after plug connect power, the indication light of adapter will on.

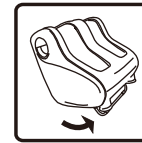

Open or close bracket according to the height of fee.

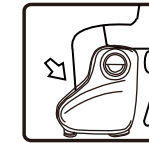

Choose chair be equal to stature, put feet into massage slot.

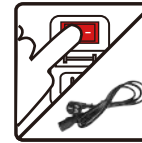

stop use: please press the off button ,tum off the power switch, then cut the power remove the adapter



If the electricity was cut off in the course of using ,step one: turn your feet side and move easily.



If the electricity was cut off in the course of using ,step one: drive up your knee.

Caution

- Do stand on , or exert excessive pressure on the massage. This may significantly slow down or stop tie motor ,potentially cause overheating and subsequent damage to the massager.
- The total recommended usage time of the massager is 15 minutes.
- Press the "stop" button to terminate massager if discomfort is experienced. The 4 pairs of discs wil not stop instantaneously ; they will continue to complete the last kneading cycle. Before stopping at their "home" position.
- Some users may experience itch on their feet during and after the massage. This a alright because it is the result of improved blood circulation.